



MICRO PULSE TECHNOLOGY

X HIFU ANTI AGING

MPT

**IMMEDIATE
RESULTS**

NO SCARS


NO DOWNTIME

**NO SIDE
EFFECTS**


SMAS FOCUS X

One Machine, Multi-Functional

For Lifting & Tightening

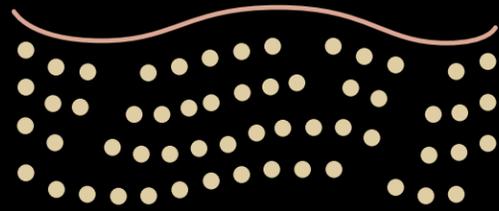


The dermis layer of the skin naturally synthesizes less collagen over the course of time.

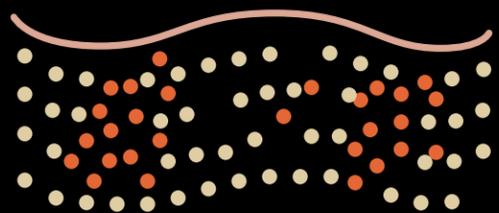

Stimulated by the focused energy, skin collagen starts to regenerate.




Focused ultrasound waves deliver energy directly within the skin layer causing a shrinking effect.

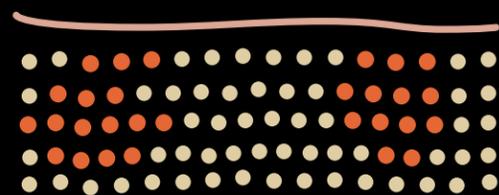


Through the regeneration of collagen, wrinkles are eliminated in the process of both firming and tightening the skin.




For Slimming & Contouring

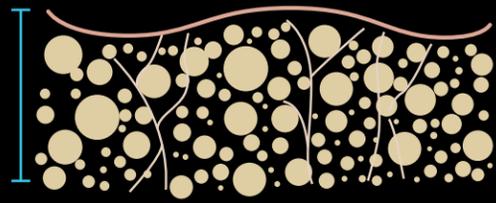

Fat cells and cellulite is formed within the body.

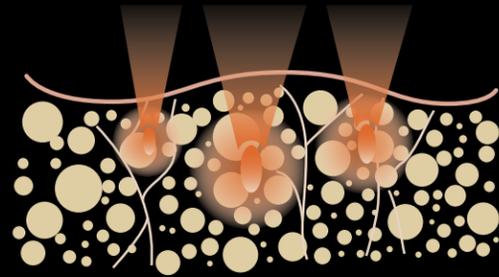


The smasfocus X energy destroys the targeted fat cells.

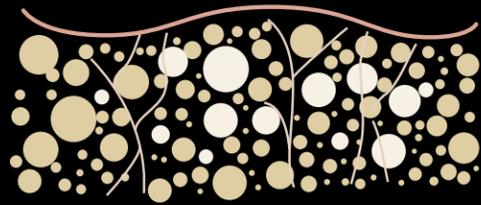

As part of the body's natural healing response, cells regenerate and join together to fill the gap.



The fat layer is naturally eliminated from the body, resulting in smooth and well-contoured body lines.

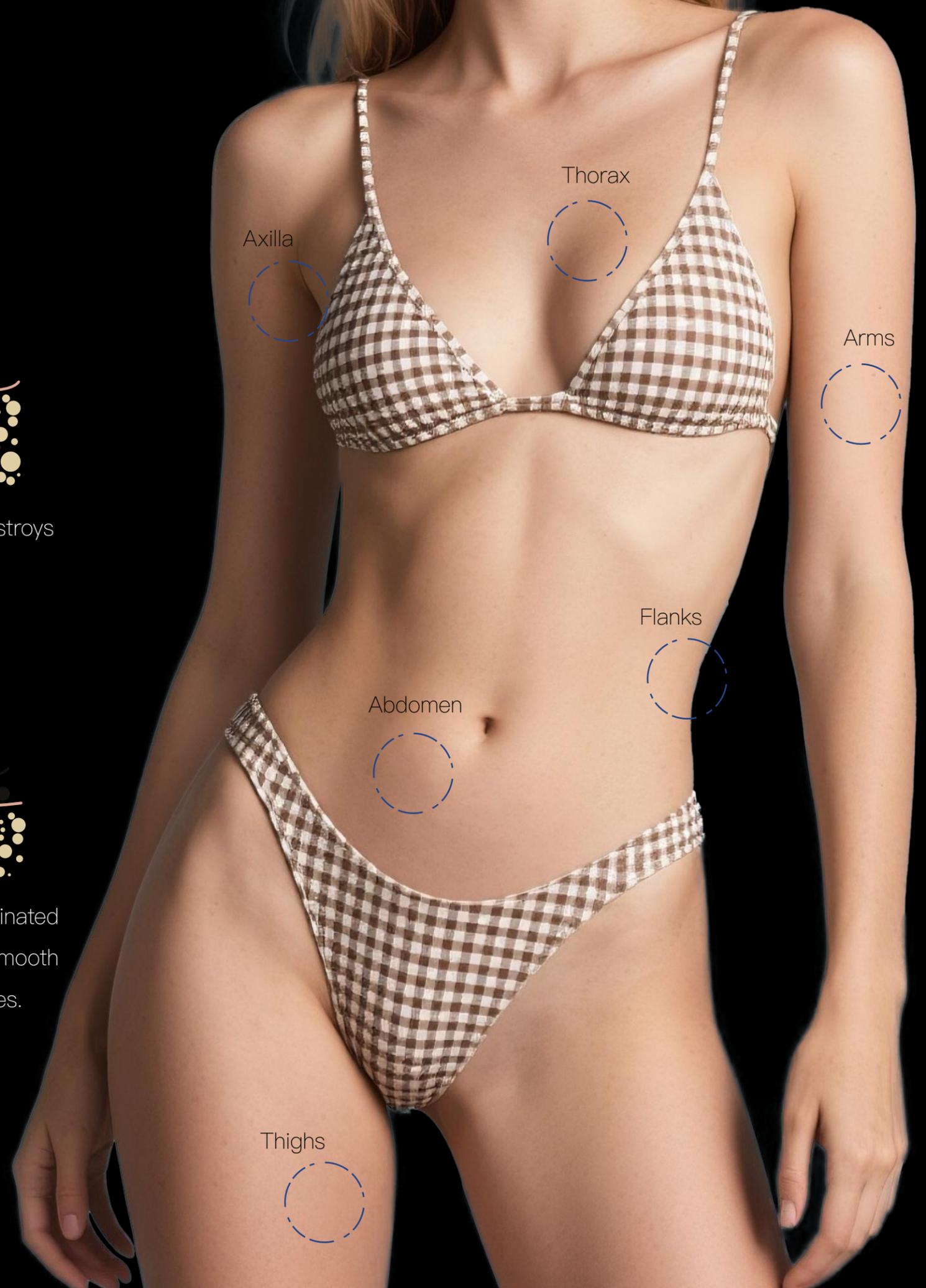

SINGLE-LINE HANDLE

7 Cartridges to Meet Various Body Part Needs

10000 Shots

FACE



1.5mm, 2.0mm, 3.0mm, 4.5mm

BODY



6.0mm, 9.0mm, 13.0mm



MICRO PULSE MODE

NORMAL MODE

Change the Mode with One Click



MULTI-LINE HANDLE

Various Cartridges Up to 8 Types

20000 SHOTS

FACE




1.5mm, 3.0mm, 4.5mm

BODY



6.0mm, 8.0mm, 10.0mm, 13.0mm, 16.0mm

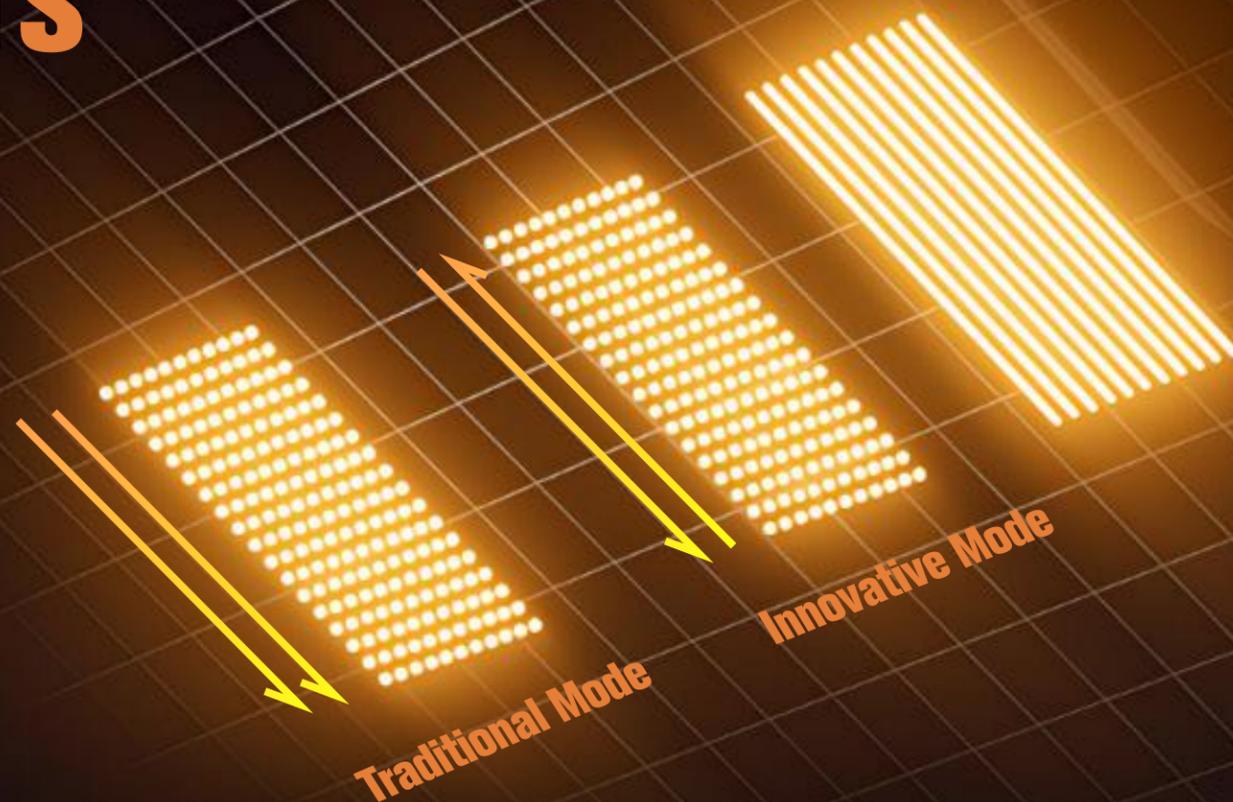

REDUCED TREATMENT TIME

The groundbreaking S-shaped working mode can reduce nearly twice the treatment time.

BOOSTER HANDLE

Frequency: 1-8hz

Making the handpiece smaller in size

Generating energy in a circular motion (90°-360° Adjustable)

60000 Shots


FACE



1.5mm,

3.0mm,

4.5mm

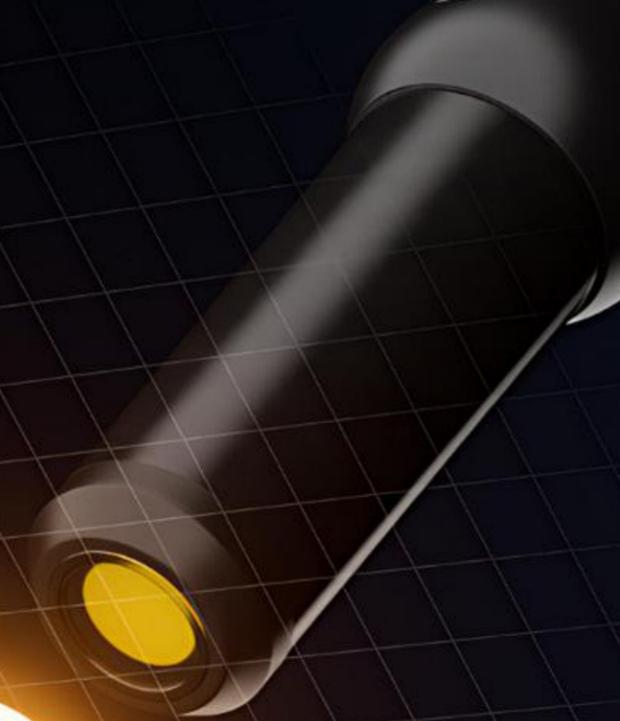

NORMAL

CIRCULAR

MICRO CIRCULAR

Double lifting effects with v_{max} overlapping treatment after the linear procedure.



Slim and light cartridges with adjustable penetration depths. Convenient and safe treatment for curved and narrow areas with better visibility on the treatment spot.

PRIVATE HANDLE (OPTIONAL)


10000 Shots



3.0mm, 4.5mm



90°-360°

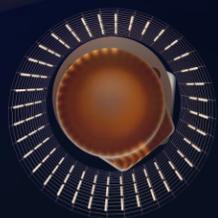

AUTOMATED OPERATION


5-90 Rotation Angle Adjustable

FOR VAGINAL TIGHTENING

Laser vs HIFU

- Collagen Rejuvenation
- Tightening

More Stronger | More Deeper
More Safer

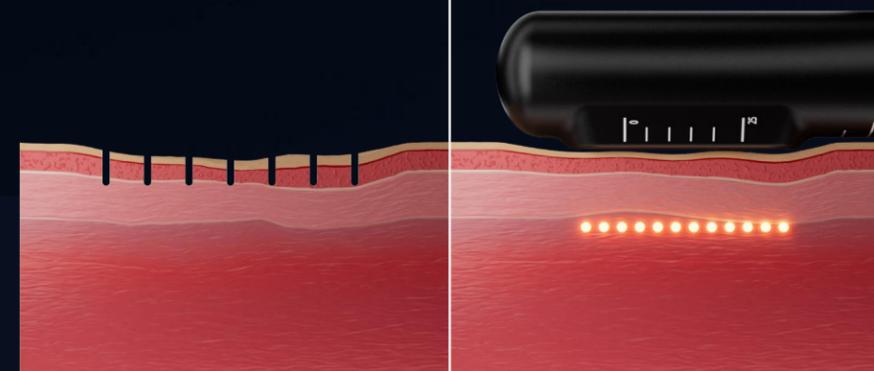

BEFORE & AFTER



Before

After



Before

After



Before

After



Before

After



Before

After



Before

After



Before

After



Before

After



Before

After


15.6-INCH HD SCREEN

Ultra-thin border ultra-large display screen ratio



Aviation-grade Power Supply

Built-in instant response sensor significantly enhances safety performance.

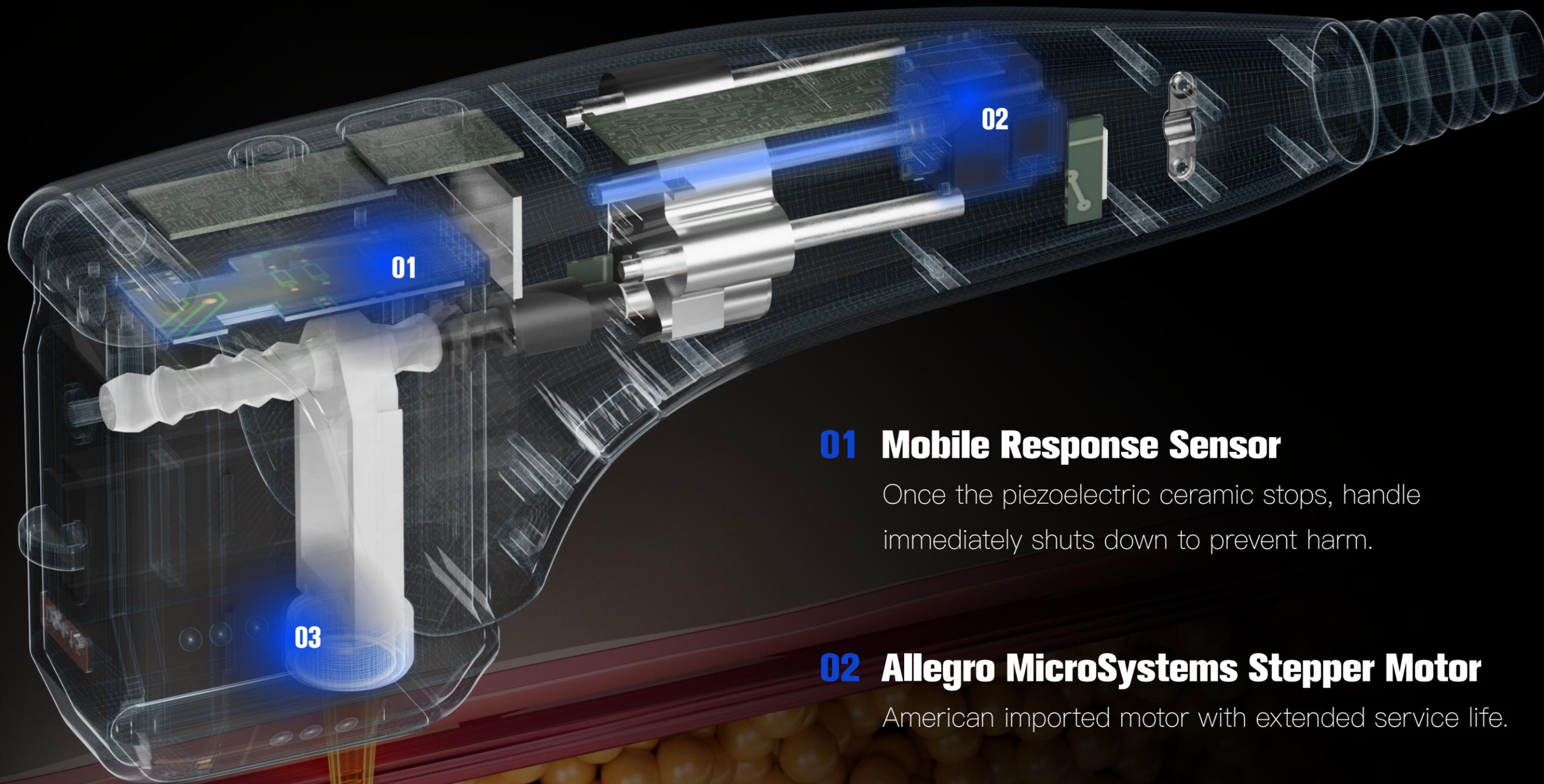

01 Mobile Response Sensor

Once the piezoelectric ceramic stops, handle immediately shuts down to prevent harm.

02 Allegro MicroSystems Stepper Motor

American imported motor with extended service life.

03 Japanese Piezoelectric Ceramics

Ensures continuous and stable energy output.

SPECIFICATIONS

Power	1200 W
Power Frequency	4mhz–7mhz
Screen	15.6 Inches
Width	1–2.5mm Adjustable
Length	8–25mm Adjustable
Dot Pitch	1.0–2.5mm Adjustable
Energy Output	0.1–2J/cm2

